


Bettersizer 2600

Particle Analysis for Your Every Need

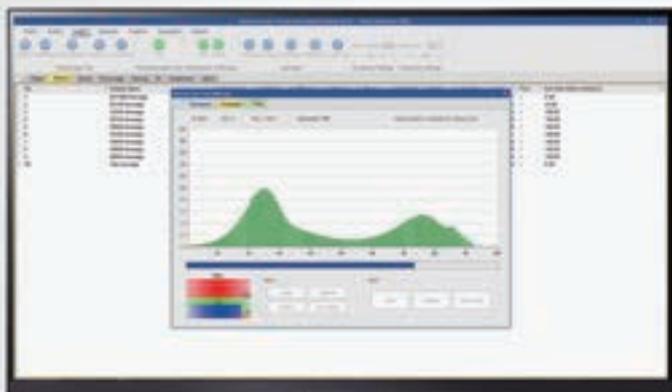

PARTICLE SIZE

PARTICLE SHAPE



Particle Analysis for Your Every Need



The Bettersizer 2600 employs a well-established particle measurement technology—Laser Diffraction Technology, covering a **particle size** range of **0.02 - 2,600 µm**. There are 92 photoelectric detectors to convert light signals from the scattering spectrum to electrical signals, which are transmitted into an intelligent software. By implementing the Mie scattering theory to the data and performing mathematical conversion, the accurate and reliable particle size distribution can be derived.

The modular design of the Bettersizer 2600 endows it with versatile capabilities. The dynamic imaging module enables combined laser and imaging tests, extending the measurement range up to **3,500 µm** and providing both **particle size and shape** analysis. The diverse dispersion system modules support a wide range of testing needs, including both dry and wet dispersion methods to meet various testing requirements.

"Ultra-Met has been using the Bettersize Laser Particle Size analyzer for over 1 year. We are able to obtain accurate, repeatable results of our tungsten carbide spray dry powder. This is a crucial tool for our company for development and quality control for our powder metal product. Customer service is always very reliable. I can always rely on having my questions answered same day by a knowledgeable member of their staff. Great company very pleased with their product and customer service overall."

Cara D'Angelo

Metallurgical Engineer
Ultra-Met Carbide Technologies, USA

Applications and Industries

| Industries | Samples | Significance |
|---|---|---|
| Pharmaceuticals  | Lactose, powder inhalers, magnesium stearate, microcrystalline cellulose, etc. | The particle size and size distribution of pharmaceuticals could affect the dissolution, body absorption, efficacy, and safety of drugs. The Bettersizer 2600 is capable of closely monitoring the particle sizes during the processes of pharmaceuticals development. |
| Abrasives  | Silicon carbide, diamond, corundum, garnet, boron nitride, etc. | Optimizing particle size and size distribution in abrasive materials maximize final products' performances while minimizing material waste. Monitoring particle size and identifying agglomerations in raw material could be achieved via the dynamic image analysis technology in Bettersizer 2600. |
| Batteries  | Lithium iron phosphate, lithium cobalt oxide, lithium manganate, modified graphite powder, etc. | The Bettersizer 2600 monitors the particle size distributions of the lithium-ion battery materials, which are critical in affecting the performances of a battery, including energy storage, stability, and safety. It is essential to strictly control the particle size distributions of lithium-ion battery materials. |
| Building Materials  | Cement, rock, clay, sand, wood, gravel, synthetic polymers, etc. | The particle size distribution of cement directly affects the hardening rate, strength, and fluidity of the final set concrete, which is the primary application of cement. Accurate and repeatable measurements of the cement particle size provided by the Bettersizer 2600 reduce costs and provide an optimized distribution in the concrete manufacturing process. |
| Paints, Inks & Coatings  | Titanium dioxide, organic pigments, iron oxides, ceramic inks, etc. | The Bettersizer 2600 characterizes particle size and size distribution of pigment-based inks. It is a crucial process in ensuring the ink remains stable over long periods of storage to prevent aggregation, color inconsistencies, and blockages in the channels or nozzles. |

| Industries | Samples | Significance |
|---|---|--|
| Mining & Minerals  | Calcium carbonate, kaolin, talc, quartz powder, graphite, barite, wollastonite, hydromagnesite, diatomite, mica, zirconium silicate, etc. | Minerals are used in many industries, including construction, fracking and abrasives. The performance of these minerals strongly depends on the size and distribution of the particles, which can be accurately measured and characterized by the Bettersizer 2600. |
| Food & Beverages  | Sugar, milk, chocolate, coffee, mayonnaise, flour, etc. | Many important characteristics of food products, such as taste, dissolution, and extraction behavior, are affected by the particle size and shape of particulate ingredients. Equipped with dynamic image analysis technology, the Bettersizer 2600 is an ideal particle analyzer for the food and beverage industry. |
| 3D Printing Materials  | Polylactic acid (PLA), acrylonitrile butadiene styrene (ABS), alloys, intermetallic compounds, ceramics, organic materials, etc. | The particle size distribution and particle shape of 3D-printing raw materials determine the degree of surface smoothness of the final printed product. Using the Bettersizer 2600 with the PIC-1, the particle size distribution and particle shape could be optimized, thereby controlling the quality and spreadability of the powders for additive manufacturing. |
| Ceramics  | Aluminum oxide, zirconium oxide, iron oxide, etc. | Particle size analysis by the Bettersizer 2600 can help the manufacturers to determine the optimum time and temperature required for the green body, as a ceramic powder with a proportion of smaller particles reduces the sintering time due to its larger surface area. The dual camera optical system in PIC-1 effectively captures the images of agglomerated oversized particles during the R&D process. |
| Cosmetics  | Lipstick, mascara, eye shadow, moisturizer, skin cream, etc. | For cosmetics, the Bettersizer 2600 aids in monitoring subtle differences in color and shine controlled by differences in the particle size distribution. The smoothness or UV light-blocking properties of creams also vary depending on the particle size distribution. |

Why the Bettersizer 2600?

The Bettersizer 2600 excels in particle analysis through its dual optical systems: laser diffraction and dynamic imaging. The combination of dual optical systems allows for comprehensive particle characterization, making the Bettersizer 2600 a versatile and indispensable tool for advanced particle analysis.

The laser diffraction system, supported by two robust patents, ensures precise and reliable particle size results in compliance with ISO 13320, enabling a wide range of industries and applications to achieve new levels of performance. Owing to its modular design, the Bettersizer 2600 effortlessly integrates a dynamic imaging system, extending the measurement range and providing individual and quantitative particle shape analysis in real time with ISO 13322-2 compliance.



Laser Diffraction System

Patented Technologies Driving Instrument Excellence



I. Combination of Fourier and Inverse Fourier Design

The Bettersizer 2600 is superior in the combination of Fourier and inverse Fourier design. Its laser system structure features 92 detectors in total, including forward, lateral, and backward detectors.

Equipped with a widely distributed spherical detector array, the Bettersizer 2600 can detect light signals across a broad angular range from 0.016° to 165°, enabling precise measurement of both small and large particles.

Total internal reflection occurs when the light transitions from a denser medium (glass) to a rarer one (air) and the incidence angle exceeds the critical angle, limiting the angles at which light can escape.

The Bettersizer 2600's innovative sample cell, with its tilted design, effectively minimizes total internal reflection. This allows more light signals to reach the detectors, enhancing measurement reliability and ensuring the acquisition of more comprehensive sample information.

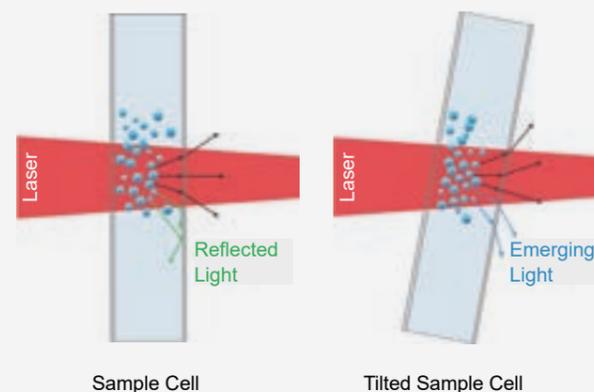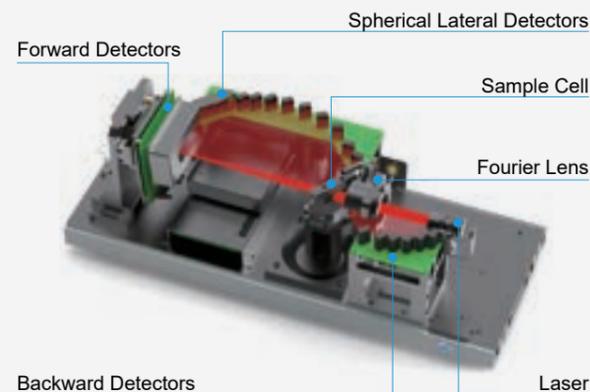

I don't know the refractive index of my sample.



No problem. With one simple click, Bettersizer 2600 can provide this parameter to you.

II. Refraction Index Measurement


Under Mie theory, measurements by laser diffraction can be particularly challenging for samples due to a variety of factors, including the following:

- Samples with completely unknown complex refractive index;
- Samples with heterogeneous chemical composition;
- Samples with significantly different particulate optical properties compared to the bulk material;
- Samples having a distinctly strong optical dispersion (small Abbe number).

To address these challenges, determining the refractive index is one of the most effective solutions. The Bettersizer 2600 offers the following capabilities:

- Determine refractive index for samples with unknown refractive index;
- Measure samples with unknown properties;
- Verify the known data of a material at a specific light wavelength;
- Provide key parameters to calculate particle size distribution in real time.

| Material | Refractive index (literature) | Refractive index (measured) |
|---------------------------|-------------------------------|-----------------------------|
| CaCO ₃ | (1.53 - 1.65) - 0.1i | 1.62 - 0.1i |
| BaSO ₄ | 1.65 - 0.1i | 1.68 - 0.1i |
| ZnO | 2.008 - 0.1i | 2.02 - 0.1i |
| Carbon black | 1.88 - 0.55i | 2.02 - 0.1i |
| Al Powder | 1.4 - 3.9i | 1.42 - 3.0i |
| SiO ₂ - Quartz | 1.54 - 0.00i | 1.54 - 0.01i |

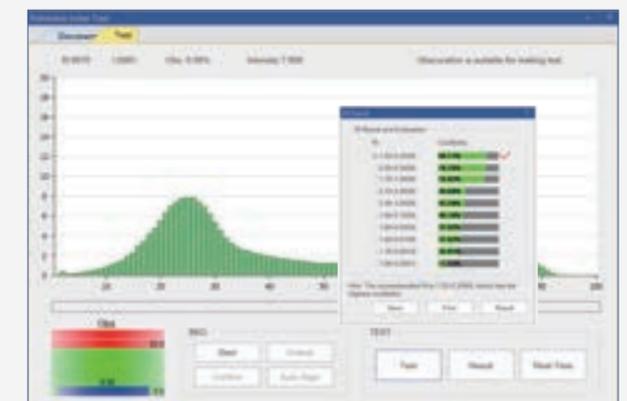

Laser Diffraction System

Superior Performance in Particle Size Analysis


I. Wide Measurement Range

Due to the instrument's excellent laser system design with 92 detectors and a very wide angular range from 0.016° to 165°, it achieves a measurement range from 0.02 µm to 2,600 µm, covering both nano and millimeter scales.




II. High Resolution

The Bettersizer 2600 is capable of distinguishing different samples with varying particle sizes within a single measurement due to its high-resolution analysis ability.

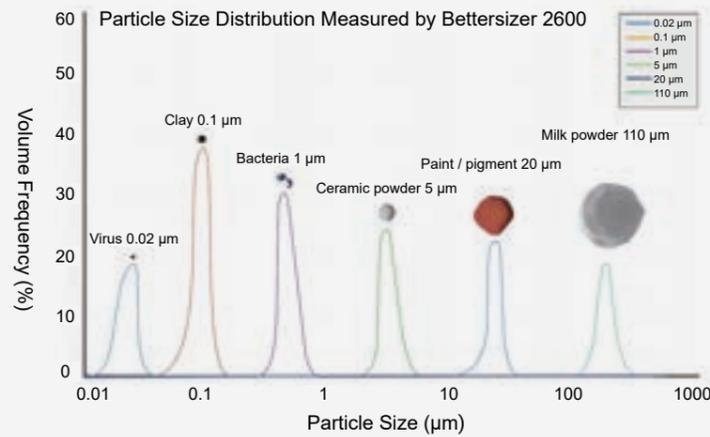

III. High Sensitivity


When gradually adding one sample to the other, the Bettersizer 2600 displays the change of particle size distributions in the curve, verifying its excellent sensitivity.

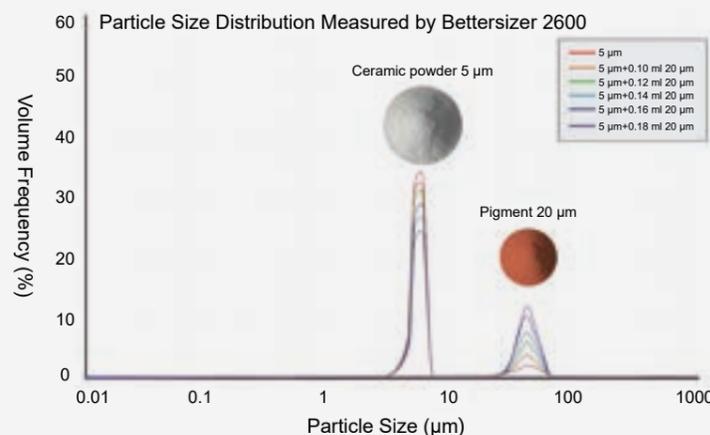

Dynamic Imaging System

PIC-1: Modular Dual-Camera Imaging System Design

Dual-Camera Imaging System

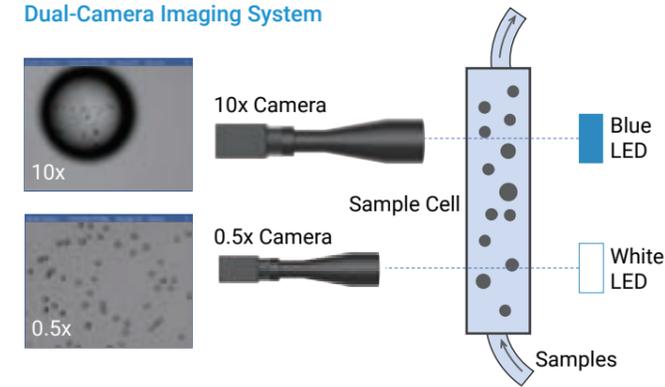


The PIC-1 features a dual-camera dynamic imaging system, making it a leading dynamic imaging module that seamlessly integrates with the Bettersizer 2600. As the dispersion system transports particles through the sample cell, the high-speed cameras capture and convert images to digital format for real-time analysis.

Going beyond mere particle size distribution, this capability allows scientists, researchers, and engineers to utilize particle shape characteristics for a deeper understanding of particles.

Broadened Insight with Modular Imaging Expertise

Extended Measurement Range

The seamless combination of laser diffraction and image analysis broadens the measurement range of Bettersizer 2600 to an impressive 3,500 µm. Samples with extremely broad distributions are now possible to measure, such as river sediment.

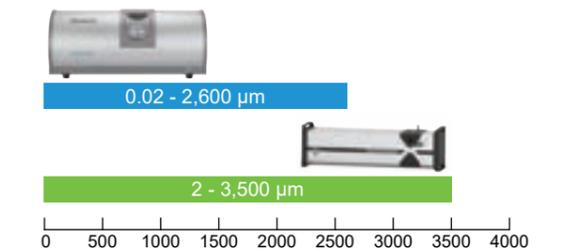


24 Detailed Particle Parameters

The PIC-1 offers comprehensive data by capturing real-time particle images, allowing customers to study individual particles in detail. With the ability to analyze 24 particle size and shape parameters, the PIC-1 provides an in-depth understanding of particle characteristics.



Oversized Particle Detection for Powder Consistency

The combination of laser diffraction and image analysis can sensitively detect oversized particles that are statistically underrepresented within a wide-distributed sample, such as oversized grain, agglomerates, air bubbles, etc.

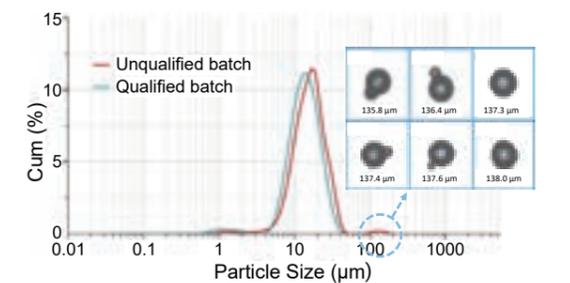

Bettersizer 2600 Family

Wet Dispersion Modules

| Parameter | BT-804 | BT-80N | BT-80N Pro | BT-802 |
|----------------|----------------|----------------|-----------------|-----------------|
| Liquid volume | ≤ 8 mL | 50 - 80 mL | 80 - 200 mL | ≤ 600 mL |
| Automation | Semi-automated | Semi-automated | Fully automated | Fully automated |
| Anti-corrosive | Yes | Yes | Yes | No |



Dry Dispersion Modules


| Parameter | BT-902 | BT-903 |
|-------------|-----------------|-----------------|
| Powder mass | 0.2 - 10 g | 0.02 - 1 g |
| Automation | Fully automated | Fully automated |

Imaging Module

| Parameter | PIC-1 |
|-------------|-------------------------|
| Principle | Dynamic image analysis |
| Application | Particle shape analysis |


Wet Dispersion System

The wet dispersion process begins with adding water or an organic solvent as the dispersing medium. A stirrer prevents settling, and an ultrasonic transducer aids dispersion and removes bubbles. The dispersed particles then flow into the optical systems for measurement and are recirculated to ensure continuous analysis.

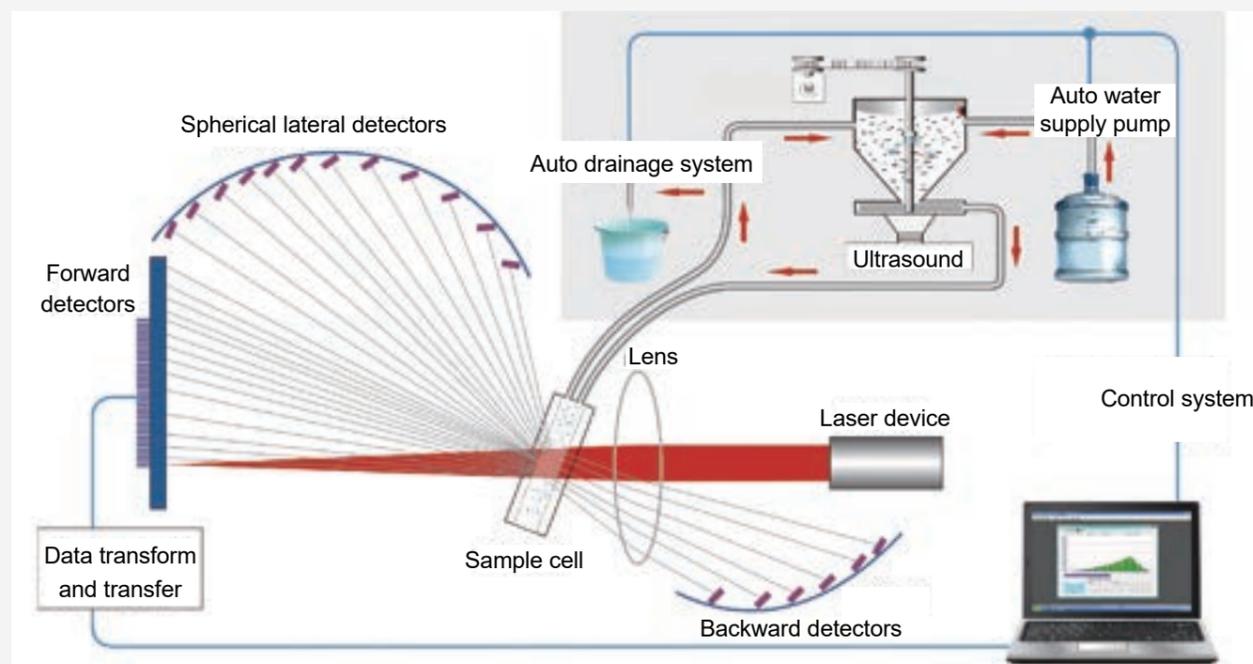

Dry Dispersion System

Dry powder samples are loaded into a sample container, which feeds them into the main disperser through vibration or high-pressure air produced by an air compressor. Inside the disperser, the particles are separated and any clumps are broken up. The dispersed particles then pass through the laser diffraction system for measurement before being collected by a vacuum cleaner.

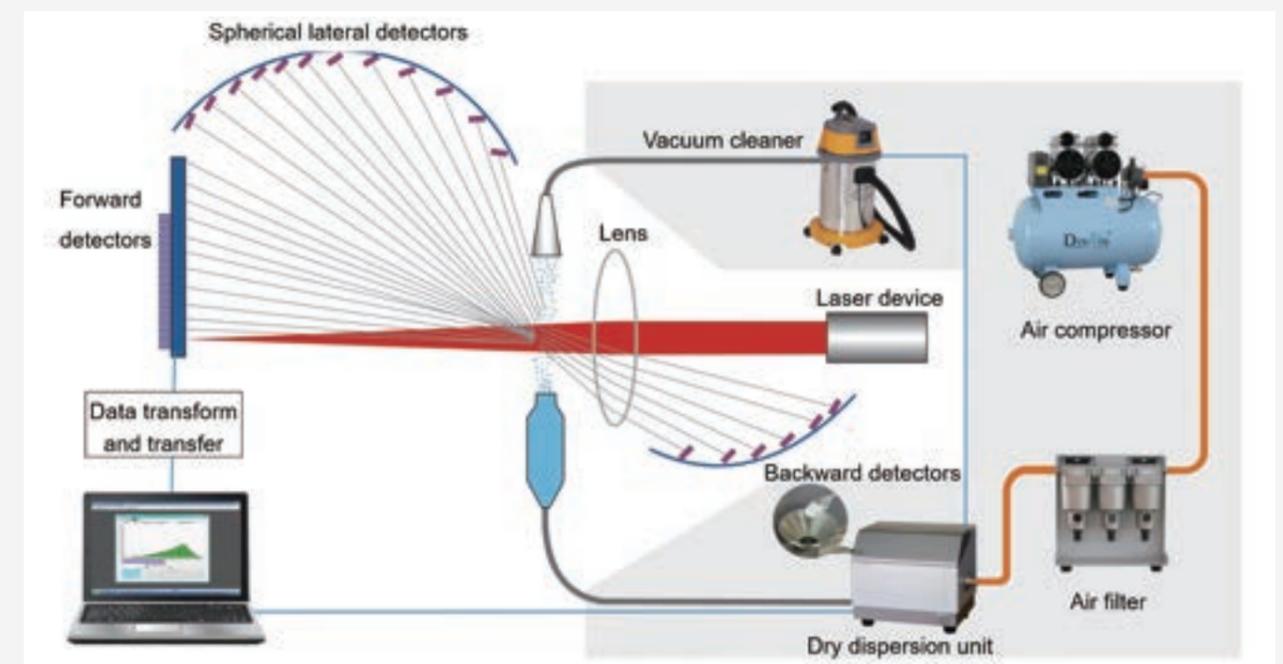

Wet Dispersion Modules

BT-802 - Automatic Wet Dispersion Unit



BT-802 is designed for particle dispersion with water as the medium. It is made up of ABS shells. The components of it include centrifugal pump, peristaltic pump, ultrasonic disperser, pinch valve, control circuit, etc.

| Parameter | Specification |
|-----------------------|--------------------|
| Measurement range | 0.02 - 2,600 µm |
| Stirring speed | 300 - 2,500 rpm |
| Ultrasonic power | 50 W max |
| Volume | 600 mL max |
| Medium | Water |
| SOP | Yes |
| Dimension (L × W × H) | 255 × 245 × 310 mm |

BT-80N and BT-80N Pro - Anti-corrosive Wet Dispersion Unit



BT-80N



BT-80N Pro

BT-80N Pro is an automatic particle dispersion unit with organic solvents. The manual BT-80N is a basic entry-level model. Both models include stainless steel shells, centrifugal pump, ultrasonic disperser, PTFE pipeline, etc.

| Parameter | Specification | |
|-----------------------|---|--------------------|
| | BT-80N | BT-80N Pro |
| Measurement range | 0.02 - 2,600 µm | |
| Stirring speed | 300 - 3,000 rpm | 300 - 2,500 rpm |
| Ultrasonic power | 50 W max | |
| Volume | 50 - 80 mL | 80 - 200 mL |
| Medium | Ethanol, Methanol, Isopropanol, Ether, Toluene, Xylene, Acetone, Octane, NMP solvents, etc. | |
| SOP | No | Yes |
| Dimension (L × W × H) | 290 × 210 × 375 mm | 240 × 220 × 290 mm |

BT-804 - Small Volume Wet Dispersion Module



BT-804 is designed for valuable or small-volume sample measurements, where the medium is water or organic solvent.

| Parameter | Specification |
|-----------------------|--------------------------|
| Measurement range | 0.02 - 2,600 µm |
| Stirring | Semi-automated |
| Volume | 8 mL max |
| Medium | Water or organic solvent |
| SOP | No |
| Dimension (L × W × H) | 200 × 160 × 200 mm |

Dry Dispersion Modules

BT-902 - Automatic Dry Dispersion Unit



BT-902 is suitable for the dispersion of dry powder particles with compressed gas. BT-902 is made up of electromagnetic vibration feeder, venturi pipe, gas circuit, electric circuit, pressure sensor, etc.

| Parameter | Specification |
|-----------------------|------------------------------|
| Measurement range | 0.1 - 2,600 µm |
| Powder mass | 0.2 - 10 g |
| Air pressure | 0.1 - 0.8 MPa |
| Funnel height | 0.7 - 2.9 mm |
| Medium | Air, nitrogen or noble gases |
| SOP | Yes |
| Dimension (L × W × H) | 310 × 190 × 255 mm |

BT-903 - Small Volume Dry Dispersion Unit



BT-903 is designed for the dispersion of small amount dry powders with a minimum sample volume of 20 mg, using compressed gas. BT-903 is composed of venturi pipe, gas circuit, electric circuit, sample tube, etc.


| Parameter | Specification |
|-----------------------|------------------------------|
| Measurement range | 0.1 - 2,600 µm |
| Powder mass | 0.02 - 1 g |
| Volume | 5 mL max |
| Air pressure | 0.1 - 0.8 MPa |
| Medium | Air, nitrogen or noble gases |
| SOP | Yes |
| Dimension (L × W × H) | 195 × 260 × 245 mm |

Effective Particle Dispersion

The Bettersizer 2600's advanced wet and dry dispersion systems are engineered to optimize dispersion energy and ensure effective dispersion without fragmentation.

Wet Dispersion: Precisely controls and monitors ultrasonication and mechanical stirring in liquid media.

Dry Dispersion: Provides exact air pressure for effective dispersion of dry powders.

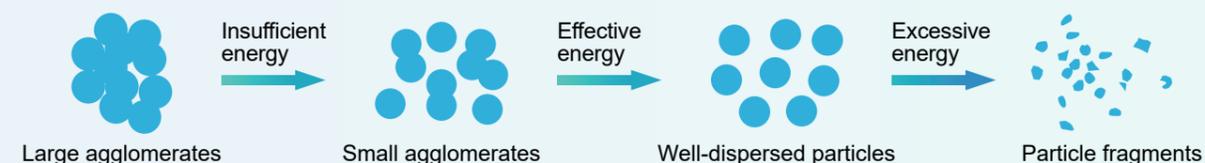

Imaging Module

PIC-1 Dynamic Imaging Module





The PIC-1 dynamic imaging module is a versatile and comprehensive dynamic image analysis instrument, designed for seamless integration with the Bettersizer 2600 laser particle size analyzer and the wet dispersion system. It primarily comprises two high-speed cameras, white and blue LED lights, and a sample cell, among other components. Utilizing dynamic image analysis, the PIC-1 precisely captures high-resolution images of particles in real time as they flow through the sample cell, allowing for detailed analysis of both particle size distribution and shape characteristics.



| Parameter | Specification |
|----------------------------------|-----------------------------------|
| Particle size range | 2 – 3,500 μm* |
| Size and shape parameters | 24 |
| Camera type | CMOS sensor, 1.5 Megapixels |
| Magnification | 0.5x and 10x |
| Imaging rate | 120 fps |
| Illumination | White and blue LEDs |
| Image recognition | Up to 10,000 particles per minute |
| Compliance | ISO 13322-2 |
| Number of size and shape classes | ≤100 (adjustable) |
| Voltage | AC 100 – 240 V, 50/60 Hz |
| Dimensions (L × W × H) | 460 × 125 × 135 mm |
| Computer configuration | |
| Processor | Intel Core i3 or a higher |
| Memory | 8 GB or higher |
| Hard disk space | 1 T or higher |
| Motherboard | with PCI-E X16 interface |
| Computer system | Windows 7 or higher |

*Wet method only



Typical Size Parameters

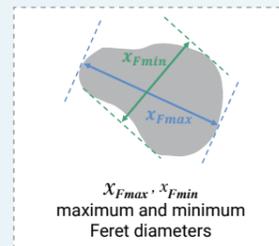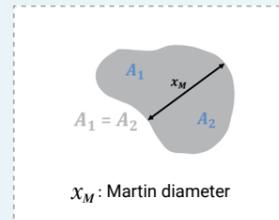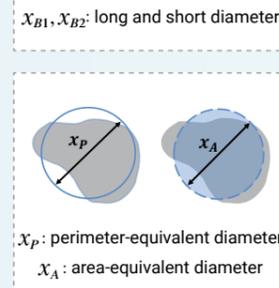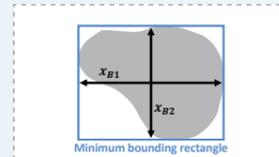

Typical Shape Parameters


$$\text{Aspect ratio} = \frac{x_{Fmin}}{x_{Fmax}}$$

$$L/D \text{ ratio} = \frac{x_{B1}}{x_{B2}}$$

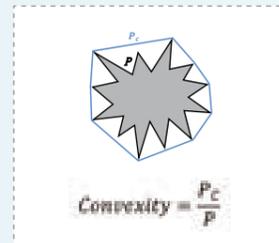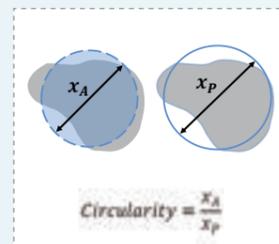

Easy Module Selection for Bettersizer 2600

The Bettersizer 2600 features both wet and dry dispersion modules for effective particle dispersion, along with a dynamic imaging module to expand the measurement range and perform particle shape analysis. Our decision tree* for module selection helps users choose the right dispersion module and determine if the dynamic imaging module is necessary.

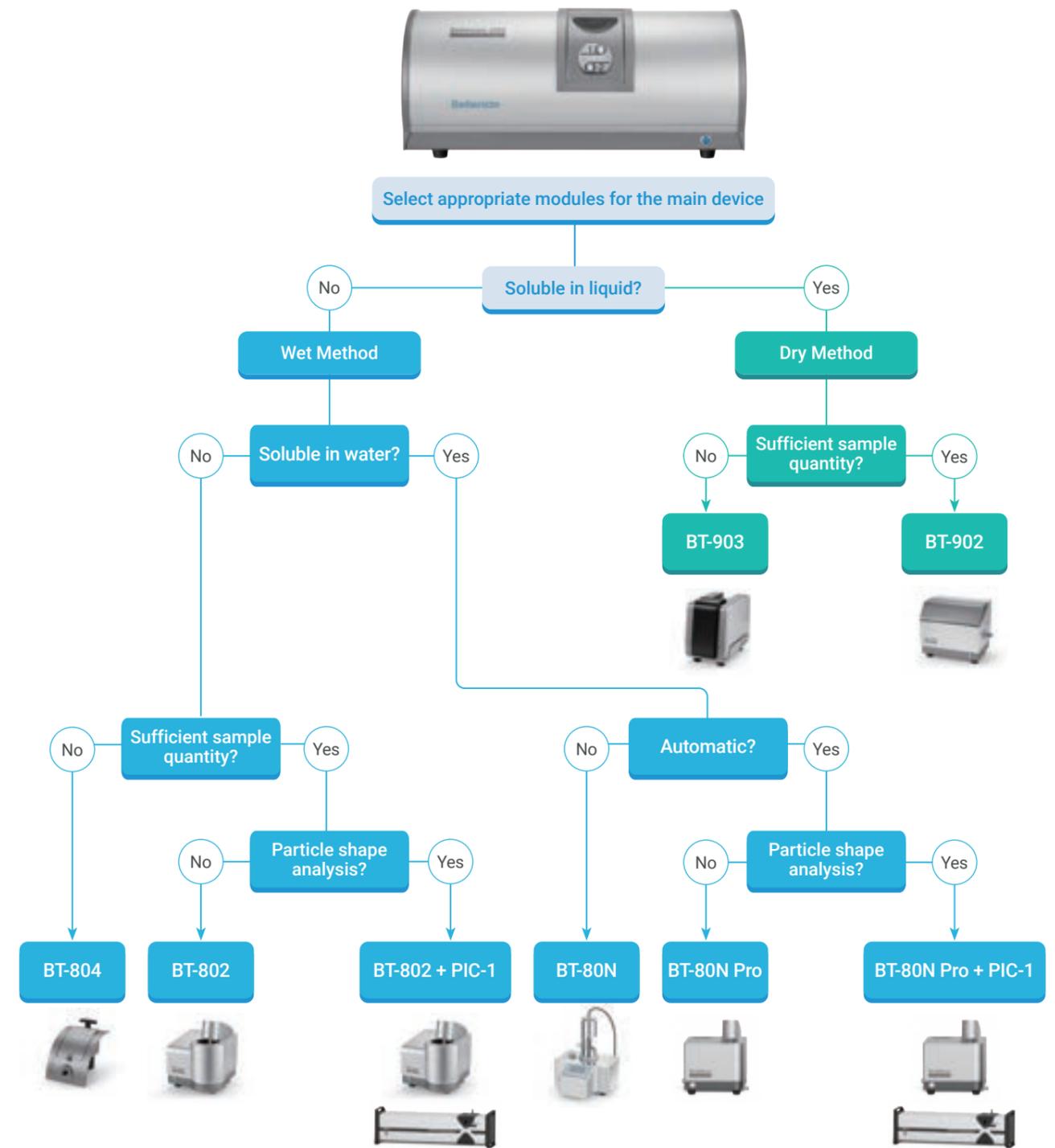

* This decision tree outlines a basic workflow for module selection. For more specialized requirements, please contact Bettersize for customized solutions.


Smart and Powerful Software

Bettersizer software is engineered to enhance your entire measurement process, from pre-processing to final data analysis. By seamlessly automating routine tasks and integrating with the Bettersizer 2600, the software empowers you to focus on scientific intricacies. Experience a streamlined workflow that delivers precise and reliable results with Bettersizer software, accelerating your research, and driving innovation.


I. Versatile Test Setup

Users can easily create new tests based on laser diffraction and dynamic image analysis methods. The software supports both automated and manual operations, providing flexibility for various sample types and testing conditions. The Standard Operating Procedure (SOP) offers a streamlined solution for standardized and automatic testing, ensuring operator-independent results that are objective and reliable.

Combined test

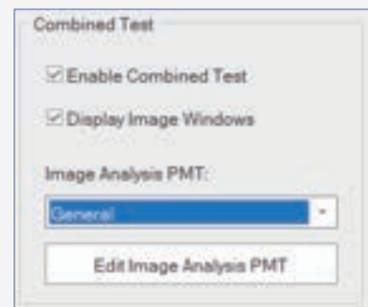


New SOP screen

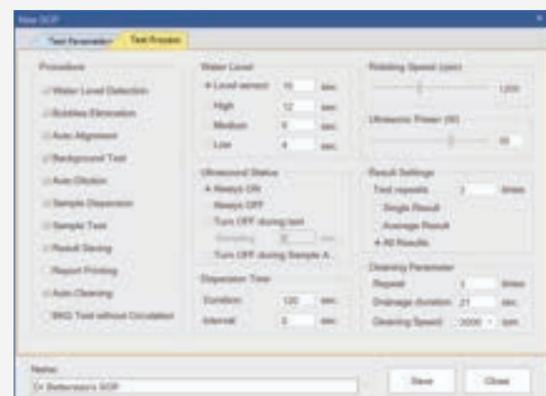


II. Automatic Pre-processing

The Bettersizer software significantly enhances data quality by automating critical instrument functions like system cleaning, optical alignment, and sample dispersion. These automated processes ensure optimal instrument performance, leading to increased precision, accuracy, and reproducibility of results.

Cleaning setting

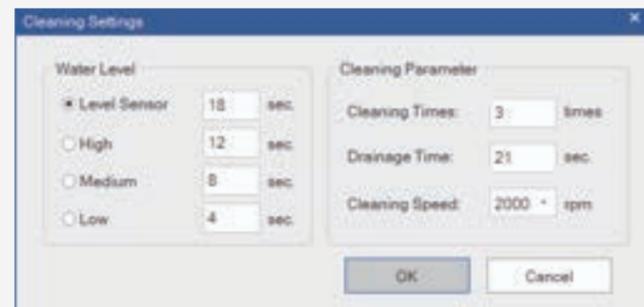


Auto alignment

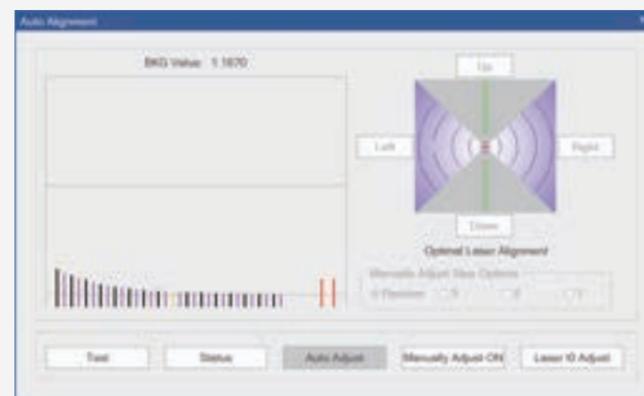

III. Real-time Testing

During the testing process, Bettersizer software delivers real-time insights into particle size distribution and shape. These immediate results provide valuable information on test progress and outcomes, enabling precise adjustments to achieve optimal results.

Imaging windows

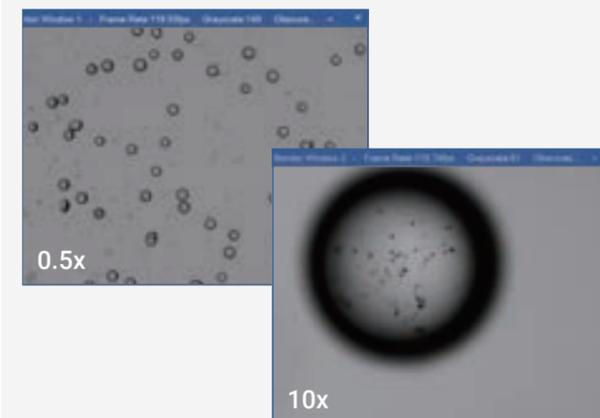

Laser diffraction test




IV. Comprehensive Data Analysis

The Bettersizer software excels in delivering comprehensive data analysis and report generation capabilities. The software allows users to customize and edit reports to meet specific requirements, including various data points, charts, and graphical representations, to create clear and informative reports. The data evaluation tools can help in assessing the result quality.

Graphical editing




Data evaluation

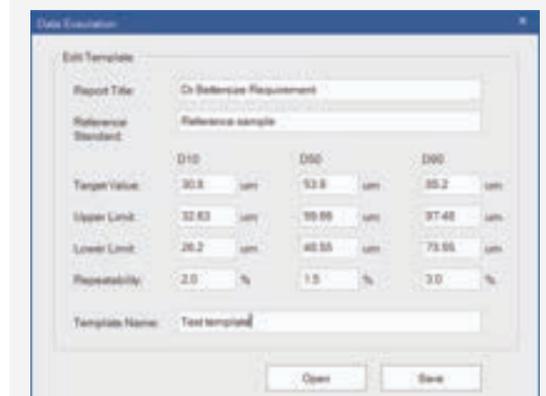

Highly Customizable Report

- **Complete and detailed data:** Frequency and cumulative distribution curves, simplified and complete distribution table, etc.
- **Editability:** Users can easily edit the reports and change the font, layout, format, etc.
- **Convertibility:** Users can switch the formats of reports freely among PDF, Excel, Text, etc.

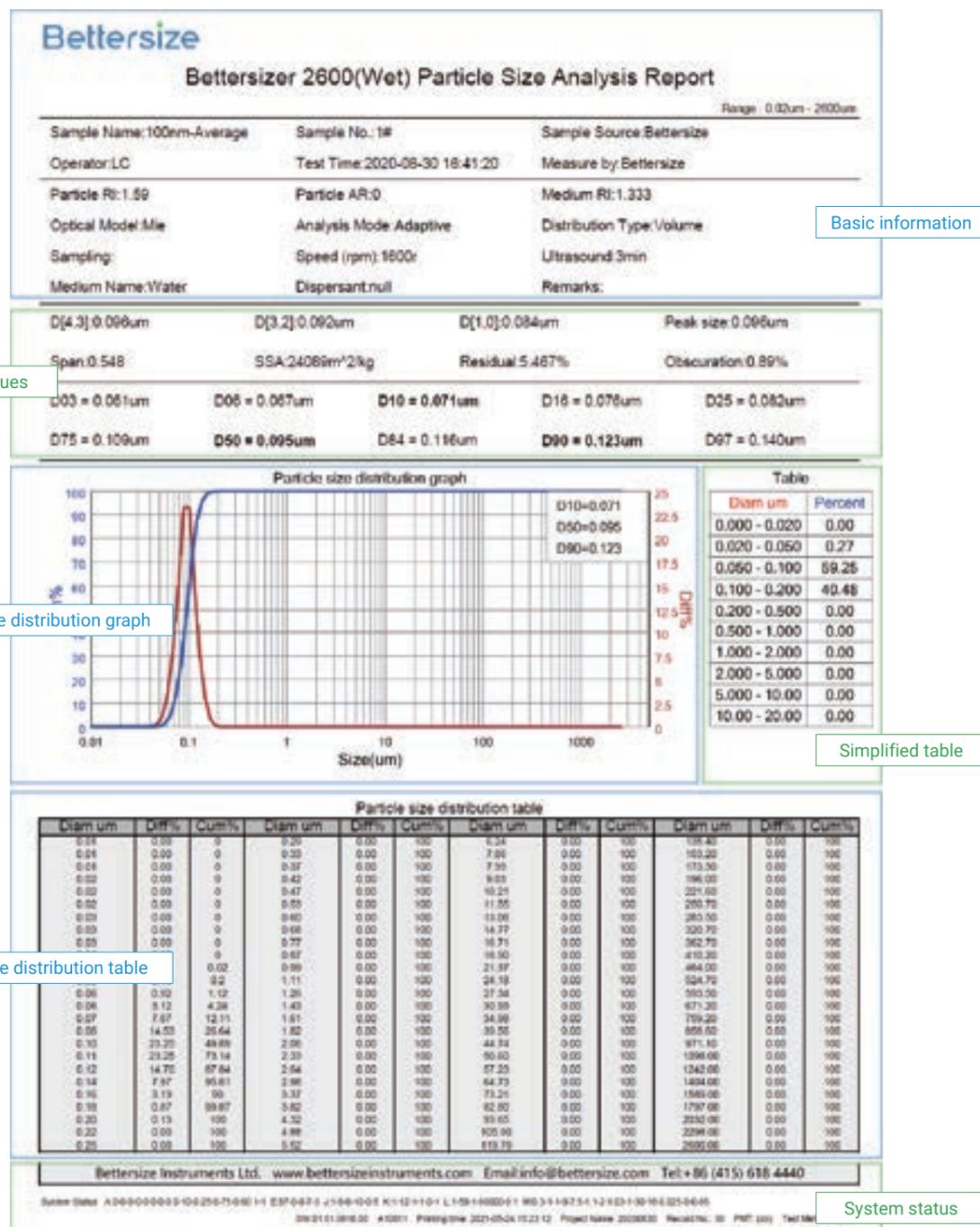

Basic information

Typical values

Particle size distribution graph

Simplified table

Particle size distribution table

System status

"Bettersizer 2600 has a superior performance on drug quality control with its rapid measurement and easy operation."

SINOWAY Pharmaceutical

Specifications

| General | |
|--------------------------------------|--|
| Principle | Laser diffraction technology; Dynamic image analysis |
| Analysis | Mie scattering theory and Fraunhofer diffraction theory; Image analysis |
| Typical measurement time | Less than 10 seconds |
| Measurement performance | |
| Measuring range | 0.02 - 2,600 μm (wet)*; 0.1 - 2,600 μm (dry)*; 2 - 3,500 μm (dynamic image)* |
| Accuracy | ≤ 0.5% * |
| Repeatability | ≤ 0.5% * |
| Number of size classes | 100 (adjustable) |
| Feeding mode | Automatic circulation or micro cuvette (wet), gas transportation (dry) |
| Special functions | SOP settings, refractive index measurement, sample ratio calculation |
| Main device | |
| Optical system | Laser diffraction system |
| Laser | 10 mW, 635 nm, Class 1 laser |
| Detector | 92 detectors |
| Measuring angle | 0.016 - 165° |
| Dynamic imaging module | |
| Optical system | Dynamic imaging system |
| CMOS camera | 0.5x and 10x |
| Measuring range | 2 - 3,500 μm |
| Frame rate | 120 fps |
| Wet dispersion module | |
| Dispersion medium | Water or organic solvents |
| Stirring speed | 300 - 2,500 rpm (BT-802, BT-80N Pro); 500 - 3,000 rpm (BT-80N) |
| Ultrasonication | Dry burning prevention, 50 W |
| Dry dispersion module | |
| Dispersion medium | Air / Nitrogen / Noble gas |
| Air pressure | 0.1 - 0.8 MPa |
| Compliance | |
| System | RoHS, CE, ISO 13320, USP <429>, ISO 13322-2 |
| Software | 21 CFR Part 11 |
| System parameters | |
| Dimensions (L x W x H) | 70.5 x 31.8 x 29.5 cm |
| Weight | 23 kg |
| Supply voltage | 100 / 240 V, 50 / 60 Hz |
| Computer configuration (recommended) | |
| Computer interface | At least one high-speed USB 2.0 or USB 3.0 port required |
| Operating system | Windows 7 or higher |
| Hardware specification | Intel Core i5 Processor, 4GB RAM, 250GB HD, Widescreen monitor |

*Sample and sample preparation dependent

Bettersize
BETTER PARTICLE SIZE SOLUTIONS



BREAKING
BOUNDARIES
SHAPING
THE FUTURE

Bettersize Instruments Ltd.

No. 9, Ganquan Road, Jinquan Industrial Park,
Dandong, Liaoning, China

Postcode: 118009

Tel: +86-755-26926582

Bettersize Inc.

3185 Airway Ave, Suite C2, Costa Mesa,
CA 92626, United States

Tel: +1 833-699-7493 (SIZE)

info@bettersize.com
www.bettersizeinstruments.com



Visit Our Bettersizer 2600 Site



Visit Our Official YouTube Channel

Disclaimer: By using or accessing the brochure, you agree with the Disclaimer without any qualification or limitation. Diligent care has been used to ensure that the information in this brochure is accurate, Bettersize Instruments Ltd. shall not be liable for errors contained herein or for damages in connection with the use of this material. The information on this brochure is presented as general information and no representation or warranty is expressly or impliedly given as to its accuracy, completeness or correctness. It does not constitute part of a legal offer or contract. Bettersize Instruments Ltd. reserves the right to modify, alter, add and delete the content outlined in the brochure without prior notice and without any subsequent liability to the company.

Copyright: © 2025 Bettersize Instruments Ltd. | All Rights Reserved
13.0109.00.03